# Alecto BC-27

## Digitale oorthermometer

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# Alecto BC-27

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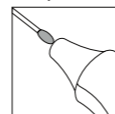
### CLEANING

#### Temperature sensor:

Clean the temperature sensor with a cotton bud soaked in medicinal alcohol (96%, available from chemists) every time after use. NEVER clean the temperature sensor with water.

#### Thermometer:

You can clean the casing of the thermometer with a damp cloth. Do not use chemical cleaning agents.

### BATTERY LOW

1. The thermometer is supplied with a 3V Lithium Batterie (CR 2032).
2. When the voltage of the battery is low, the LCD will display: . Please replace the battery in the unit. The thermometer cannot work accurately under the condition of low voltage.
3. Remove the battery cover and pull down the metal button (see Fig), the battery will be bounced out, then you can take out the old battery.
4. Place a new battery, the positive pole (+) up and push the battery down so that it clicks into place.
5. Place the battery cover back onto the thermometer.
6. Also take the battery out if the unit will not be used for over six months. This to prevent any damage in case the battery starts leaking.

### Memory:

The thermometer remembers the temperature of the last 10 measurements.

1. Switch the thermometer on by briefly pressing the **POWER/M** button.
2. Press the **POWER-M** button repeatedly to view the memories. (Memory 10 is the last measurement, Memory 09 is the last but one, etc.)
3. After 1 minute the thermometer switches off automatically.

### Possible causes of the temperature measured being too low:

- The thermometer is not far enough or at an angle in the ear. Place the thermometer in the ear correctly and measure the temperature again.
- The temperature sensor is dirty. Clean it and measure again.
- The ear canal is blocked or dirty. Clean it or measure the temperature in the other ear.

For further questions and/or problems, contact the Alecto customer service.

### Advice and warnings: Never:

- Never measure your temperature with the thermometer if you have an ear infection or other ear irritation.
- Never measure your temperature within 30 minutes of physical activity or bathing.
- Never touch the thermometer sensor.

### Always:

- Clean the thermometer well after use and store it in a dry and dust-free environment and out of reach of children.

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Display-Auflösung: 0.1°C (0.1°F)

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### ATTENTION

- Keep the thermometer and battery out of the reach of children.
- Do not talk or walk around while taking your temperature.
- Store the thermometer in the special case.
- Do not expose the thermometer to direct sunlight, dust or moisture.
- Do not dismantle the thermometer.

Ensure that button cell batteries are NEVER swallowed by babies or young children. Keep (button cell) batteries far away from babies and children at all times and always discard used batteries at a dedicated collection point. If you have any doubt about your baby or child having swallowed a button cell battery, immediately consult your doctor.

### SPECIFICATIONS

Range of displayed temperature: 34.0°C ~ 44.0°C (93.2°F ~ 111.2 °F)

Operating ambient:  
Temperature: 16°C ~ 35°C (\*) (60.8°F ~ 95.0°F)  
Humidity: 80%RH  
\*: in case the temperature is out of this range, the display shows "ER 1".

Storage and transport condition:  
Temperature: -25°C to 55°C (-13°F to 131°F)  
Humidity: 85%RH

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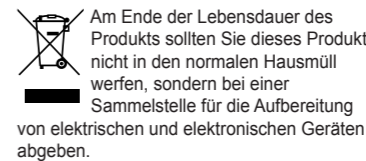
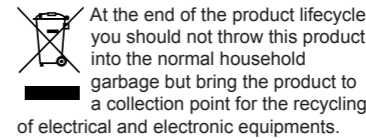
Battery:  
1x 3V Lithium battery (CR 2032)  
Battery life: 2000 measurements  
Power consumption: 9 mW

Dimension: 110mm by 35mm by 35mm  
Net weight: Approx. 37g (battery incl.)

### ENVIRONMENT

Do not throw used batteries away. Hand them in at your local depot for domestic chemical waste.

### CE

The declaration of conformity is available at:  
<http://DOC.hesdo.com/BC-27-DOC.pdf>

### WARRANTY

For the Alecto BC-27 thermometer you have a warranty of 24 months from the date of purchase. During this period, we guarantee the free repair of defects caused by material and workmanship errors. All this subjected to the final assessment of the importer.

HOW TO HANDLE: If you notice any defect, first refer to the user's manual. Please contact the supplier of this thermometer or contact customer service of Alecto via the Internet [www.alecto.nl](http://www.alecto.nl)

THE WARRANTY BECOMES NULL AND VOID: In case of improper use, incorrect connections, leaking and/or incorrectly installed batteries, use of unoriginal parts or accessories, negligence and in case of

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Hangzhou Hua'an Medical & Health Instruments Co., Ltd,  
Building 2, 1# Fuzhu Nan RD,  
Wuchang Town, Yuhang District,  
Hangzhou, Zhejiang 310023

Shanghai International Holding Corp. GmbH,  
Eiffestrasse 80,  
20537  
Hamburg Germany

Imported by: Hesdo bv.  
Azielaan 12, 5232BA  
's-Hertogenbosch  
The Netherlands

defects caused by moisture, fire, flooding, lightning and natural disasters.

In case of unauthorised modifications and/or repairs performed by third parties. In case of incorrect transportation of the device without adequate packaging and when the device is not accompanied by this warranty card and proof of purchase. All further liability, especially concerning any consequential damages, is excluded.

Hangzhou Hua'an Medical & Health Instruments Co., Ltd,  
Building 2, 1# Fuzhu Nan RD,  
Wuchang Town, Yuhang District,  
Hangzhou, Zhejiang 310023

Shanghai International Holding Corp. GmbH,  
Eiffestrasse 80,  
20537  
Hamburg Germany

Imported by: Hesdo bv.  
Azielaan 12, 5232BA  
's-Hertogenbosch  
The Netherlands

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Shanghai International Holding Corp. GmbH, Eiffestrasse 80, 20537 Hamburg Germany